

Dobot Magician E6 — quick guide

The arm is position controlled and does not feel you. It will not stop for your hand, the bench, or the phantom.

- Hold **L1** or nothing moves. Release it and the arm stops — that is your first reflex, not the emergency stop.
- Keep the **red emergency stop** on the base within reach. Press to stop, rotate to release.
- The **robot base is the centre** of the workspace: the arm reaches 450 mm in every direction, **front and back**. Nobody puts a hand inside it while it is running.

Start up

1. Power on the robot

Check the red emergency stop is released, switch on at the wall, then **tap** the power button on the base. Do not hold it.

 Wait **30-60 s** until the light is **steady blue**.

2. Connect the network

```
nmcli con up dobot-e6
```

3. Load the environment

```
cd "/home/yz22/Documents/code/Dobot arm"  
source setup_dobot.bash
```

4. Send the arm home

```
python3 src/dobot_e6_hw/scripts/go_home.py --speed 5
```

 Wait until it prints **Arm at home pose**.


5. Start the program

```
ros2 launch dobot_e6_hw real_hw.launch.py robot_ip:=192.168.5.1
```

 Wait about **15 s** until the light is **flashing green slowly**, then **25 s more** — until **Calibration complete** appears — before touching the controller. About **40 s** in total.

Stop: press **Ctrl+C** in the same terminal. The light returns to steady blue.

Joystick

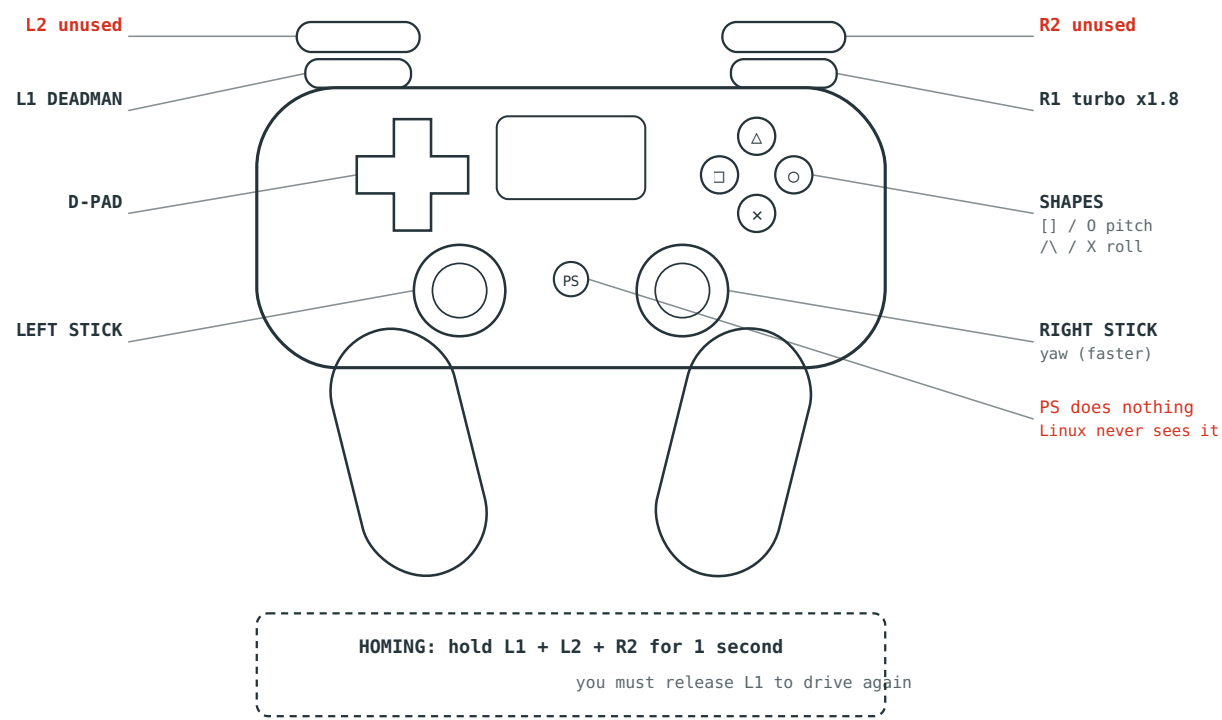

CONTROL	MOTION
L1	Hold to move. Release and the arm stops
D-pad ◀▶	Move to your left / right
D-pad ▲▼	Move away from / toward you
Left stick ▲▼	Move up / down
Right stick ▲▼	Yaw — twist
◻ / ◯	Pitch — tilt
△ / ×	Roll — rotate
L1 + L2 + R2, hold 1 s	Home — then release L1 and press it again to drive
R1	Turbo ×1.8 — operator only